


## PRODUCT DATA SHEET

### ECO ET 12 S

The LAUDA ECO stands out due to its range of functionality and ease of use. A menu navigation in plain text allows easy operation of the devices. Both control heads are equipped with a mini-USB interface as standard. Further interfaces are available as modules. Another innovation is the practical allocation of a flow-rate switch at the front on the control head. This means that individual adjustment of the flow between internal and external circulation is possible even during operation. Variants in 230 V; 50 Hz are also available with natural refrigerants. The most powerful units are equipped with the energy-saving LAUDA SmartCool system. All cooling thermostats are available as air-cooled or water-cooled variants.

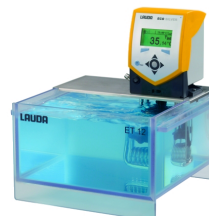

Working temperature min.  
20 °C



Working temperature max.  
100 °C

#### Technical Attributes

Technical Attributes	ECO ET 12 S
Working temperature min.	20 °C
Working temperature max.	100 °C
Ambient temperature min.	5 °C
Ambient temperature max.	40 °C
Temperature stability	0.01 ±K
Application	external   External and inside the bath   Inside the bath
Heater power	2 kW
Power consumption	2.1 kW
Bath volume min.	9.5 L
Bath volume max.	12 L
Bath opening width	300 mm
Bath opening depth	175 mm
Bath depth	160 mm
usable depth	140 mm
Dimensions (WxDxH) in mm	322x331x389

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Pfarrstraße 41/43 · 97922 Lauda-Königshofen  
Postfach 1251 · 97912 Lauda-Königshofen · DE

T +49 (0) 9343 503-0 · F +49 (0) 9343 503-222  
info@lauda.de · www.lauda.de  
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim · HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim · HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser

## PRODUCT DATA SHEET

### ECO ET 12 S

Technical Attributes	ECO ET 12 S
Pump pressure max.	0.55 bar
Pump flow max. (pressure)	22 L/min

LAUDA DR. R. WOBSE R GMBH & CO. KG  
Pfarrstraße 41/43 · 97922 Lauda-Königshofen  
Postfach 1251 · 97912 Lauda-Königshofen · DE

T +49 (0) 9343 503-0 · F +49 (0) 9343 503-222  
info@lauda.de · www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim · HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSE R Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim · HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser